


Neurosurgery

GAV[®] 2.0

Gravitational valve for the treatment
of hydrocephalus



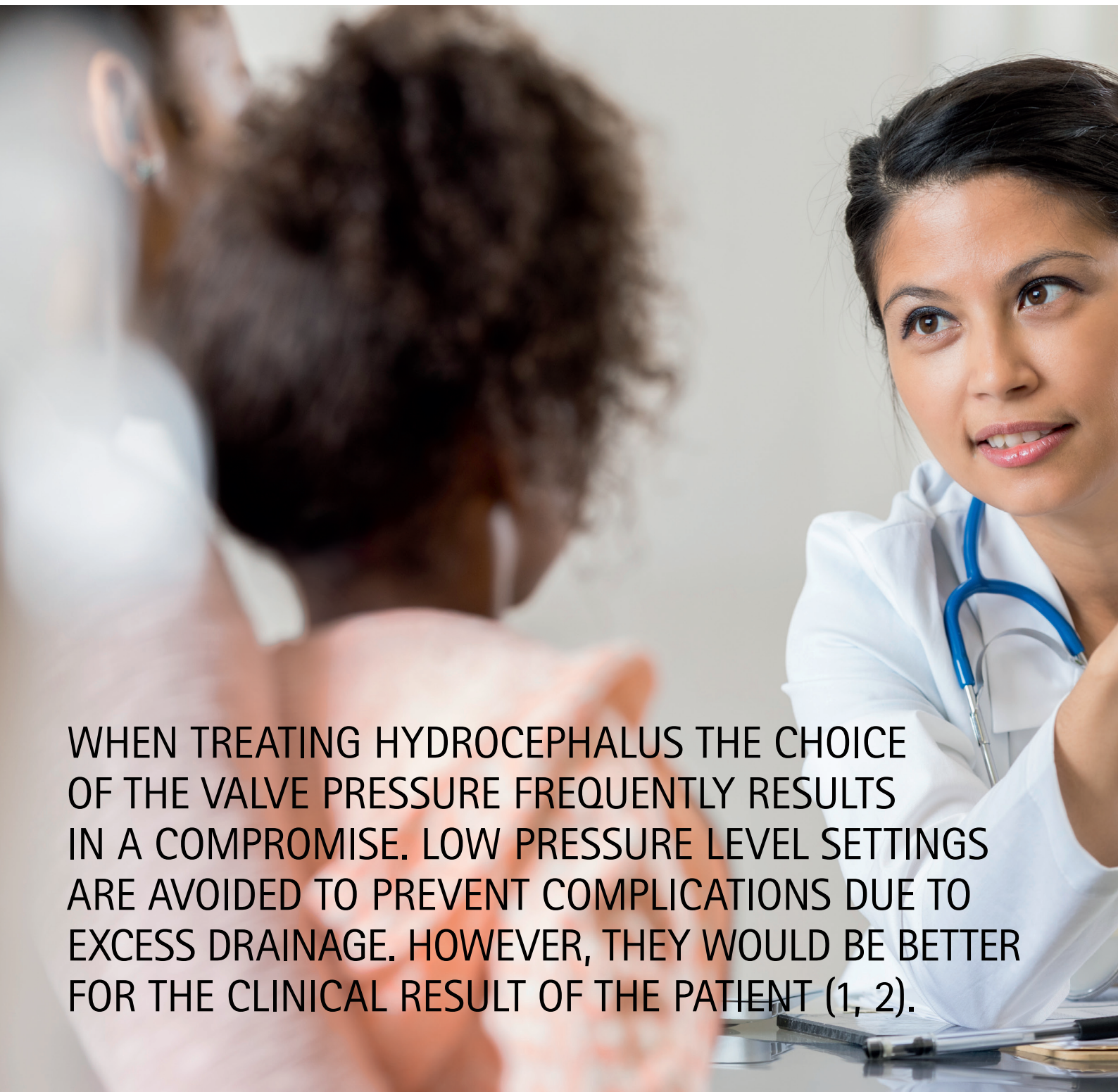

WHEN TREATING HYDROCEPHALUS THE CHOICE OF THE VALVE PRESSURE FREQUENTLY RESULTS IN A COMPROMISE. LOW PRESSURE LEVEL SETTINGS ARE AVOIDED TO PREVENT COMPLICATIONS DUE TO EXCESS DRAINAGE. HOWEVER, THEY WOULD BE BETTER FOR THE CLINICAL RESULT OF THE PATIENT (1, 2).



Conventional valve technology provides inadequate protection against overdrainage complications such as hygroma, hematoma or slit ventricles.

Excessive drainage induced by the hydrostatic suction of the vertical shunt system is considered to be one of the major causes of shunt dysfunction in the treatment of pediatric hydrocephalus (3).



Many HC-patients suffering from chronic headaches due to excessive drainage, frequently develop an irreversible slit ventricle syndrome (4, 5).

- (1) Lemcke J, Meier U, Müller C, Fritsch MJ, Kehler U, Langer N, Kiefer M, Eymann R, Schuhmann MU, Speil A, Weber F, Remenez V, Rohde V, Ludwig HC, Stengel D. Safety and efficacy of gravitational shunt valves in patients with idiopathic normal pressure hydrocephalus: a pragmatic, randomised, open label, multicentre trial (SVASONA). *J Neurol Neurosurg Psychiatry*. 2013 Aug; 84(8):850-7.
- (2) Suchorska B, Kunz M, Schniepp R, Jahn K, Goetz C, Tonn JC, Peraud A. Optimized surgical treatment for normal pressure hydrocephalus: comparison between gravitational and differential pressure valves. *Acta Neurochir (Wien)*. 2015 Apr;157(4):703-9.
- (3) Gruber RW, Roehrig B. Prevention of ventricular catheter obstruction and slit ventricle syndrome by the prophylactic use of the Integra antisiphon device in shunt therapy for pediatric hypertensive hydrocephalus: a 25-year follow-up study. *J Neurosurg Pediatr*. 2010 Jan;5(1):4-16.
- (4) Reikate HL. Shunt-related headaches: the slit ventricle syndromes. *Childs Nerv Syst*. 2008 Apr;24(4):423-30.
- (5) Buxton N, Punt J. Subtemporal decompression: the treatment of noncompliant ventricle syndrome. *Neurosurgery*. 1999 Mar;44(3): 513-8.

GAV[®] 2.0

THE VALVE

GRAVITATIONAL TECHNOLOGY

The combination of a differential pressure unit and gravitational unit ensures an automatic opening pressure adjustment as a function of the position of the patient's body and in that way counteracts complications due to excess drainage.



DESIGN

The slender cylindrical design allows for quick and easy implantation in adults and in the treatment of pediatric hydrocephalus.

GAV® 2.0 LP



GAV® 2.0 LP WITH DEFLECTION



TWO ADDITIONAL LP-VARIANTS


GAV® 2.0 is also suitable for lumbar drainage due to two special variants.

TITANIUM

The valve material titanium is biocompatible and durable. It prevents external and subcutaneous pressure influences and is MRI compatible.

GAV[®] 2.0

FUNCTIONALITY AND POSITION OF THE BODY

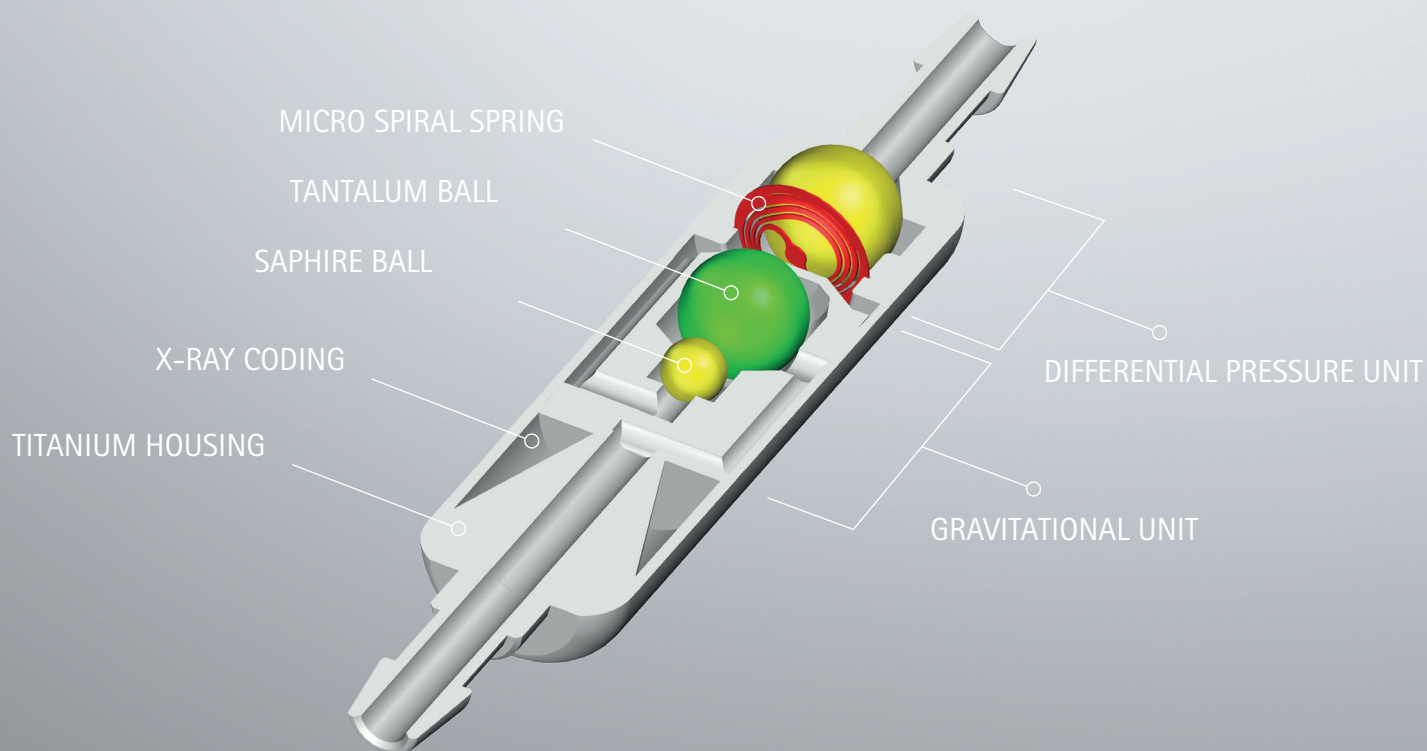

The functionality of the GAV[®] 2.0 is illustrated interactively in the Miethke App.





for Apple





for Android



Valve opening pressure
(cmH₂O)

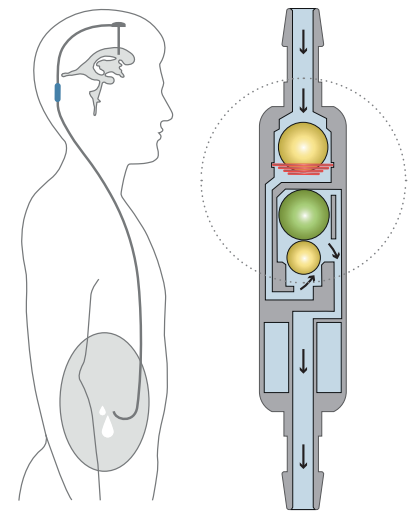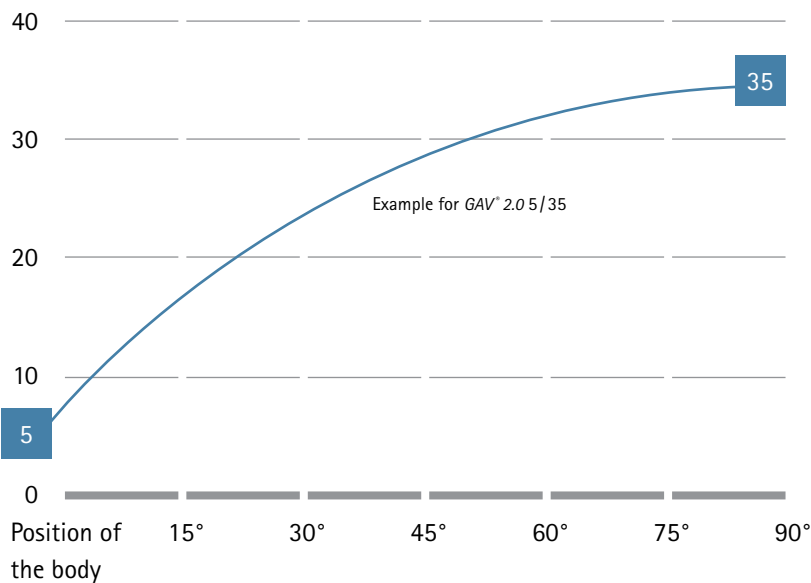

Differential pressure unit
& gravitational unit
interact in the upright
position

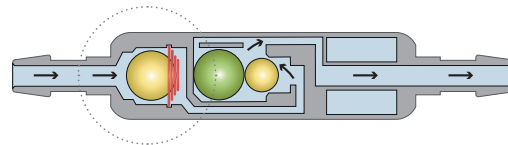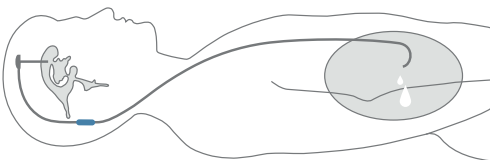

In the recumbent position, only the differential pressure unit is effective

HORIZONTAL POSITION OF THE BODY

The valve opening pressure of GAV* 2.0 in the lying position is solely determined by the micro spiral spring of the differential pressure unit. The gravitational unit is not active in this body position and is always open. If the patient's intraventricular pressure (IVP) exceeds the opening pressure of the micro spiral spring, the closure ball moves out of the cone, opening a gap for drainage purposes.

For the example, a differential pressure unit of 5 cmH₂O was selected.

VERTICAL POSITION OF THE BODY




In the vertical position of the body, the gravitational and differential pressure units act in conjunction. When the patient stands up, the tantalum ball (green) in the gravitational unit is activated, causing gravity to increase the valve opening pressure. Now, the weight of the tantalum ball (opening pressure of the gravitational unit) must be overcome in addition to the opening pressure of the differential pressure unit. Only when the sum of intraventricular pressure (IVP) and hydrostatic pressure exceeds the opening pressure of both units, drainage can take place again. The opening pressure in the upright position of the patient is therefore computed from the sum of the differential pressure and gravitational pressure.

For the example, a gravitational unit having 30 cmH₂O was selected. The total opening pressure in the upright position thus adds up to 35 cmH₂O.

GAV[®] 2.0

X-RAY RECOGNITION AND PRESSURE LEVEL RECOMMENDATION

PRESSURE LEVEL VARIANTS

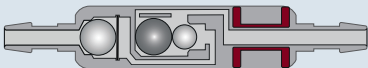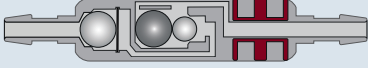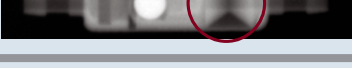Lying	Upright	X-ray coding	Radiograph
5 cmH ₂ O	20 cmH ₂ O		
5 cmH ₂ O	25 cmH ₂ O		
5 cmH ₂ O	30 cmH ₂ O		
5 cmH ₂ O	35 cmH ₂ O		
10 cmH ₂ O	25 cmH ₂ O		
10 cmH ₂ O	30 cmH ₂ O		

PRESSURE LEVEL RECOMMENDATION *

NEW BORN	CHILDREN > 6 MONTHS	CHILDREN > 3 YEARS	ADULTS	ADULTS > 65 YEARS
5/25	10/25	10/30	5/30 5/25 < 1.60 m 5/35 > 1.80 m	5/25 5/20 < 1.60 m 5/30 > 1.80 m

* Recommended pressure level in cmH₂O.

This is a non-binding recommendation. The doctor will decide each case on an individual basis.

The choice of the appropriate pressure level of GAV® 2.0 depends on several other factors, including age, degree of activity, size and stature of the patient.

The values given apply to mobile patients. For patients with little mobility or a high BMI, the gravitational unit should be chosen lower than recommended here.

GAV[®] 2.0

GAV[®] 2.0 VALVE

■ GAV[®] 2.0 Valve



13.4 mm


Valve: $d_o = 4.2$ mm
Connector: $d_o = 1.9$ mm
preferably used with
Catheter: $d_i = 1.2$ mm, $d_o = 2.5$ mm

Art. No.	Lying	Upright
FX210T	5 cmH ₂ O	20 cmH ₂ O
FX211T	5 cmH ₂ O	25 cmH ₂ O
FX212T	5 cmH ₂ O	30 cmH ₂ O
FX213T	5 cmH ₂ O	35 cmH ₂ O
FX214T	10 cmH ₂ O	25 cmH ₂ O
FX215T	10 cmH ₂ O	30 cmH ₂ O

For pressure level recommendation see page 9.

GAV® 2.0 VALVE WITH DISTAL CATHETER

- GAV® 2.0 Valve
with distal catheter (1200 mm)

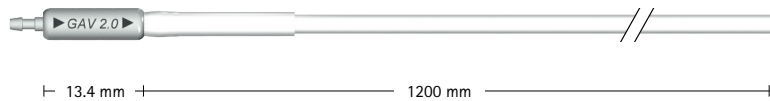

Valve: $d_o = 4.2$ mm
 Connector: $d_o = 1.9$ mm
 Catheter: $d_i = 1.2$ mm, $d_o = 2.5$ mm

Art. No.	Lying	Upright
FX216T	5 cmH ₂ O	20 cmH ₂ O
FX217T	5 cmH ₂ O	25 cmH ₂ O
FX218T	5 cmH ₂ O	30 cmH ₂ O
FX219T	5 cmH ₂ O	35 cmH ₂ O
FX220T	10 cmH ₂ O	25 cmH ₂ O
FX221T	10 cmH ₂ O	30 cmH ₂ O

For pressure level recommendation see page 9.

GAV[®] 2.0



GAV[®] 2.0 SHUNT SYSTEM WITH CONTROL RESERVOIR

- GAV[®] 2.0 Valve with integrated CONTROL RESERVOIR* and distal catheter (1200 mm)

- Ventricular catheter with pediatric deflector and introducing stylet (250 mm)

* An additional valve in the base of the CONTROL RESERVOIR makes it possible to flush the fluid only in the distal direction. This feature allows for controlling the patency of the ventricular catheter and the distal drainage.

Valve: $d_o = 4.2\text{ mm}$
Connector: $d_o = 1.9\text{ mm}$
Catheter: $d_i = 1.2\text{ mm}$, $d_o = 2.5\text{ mm}$

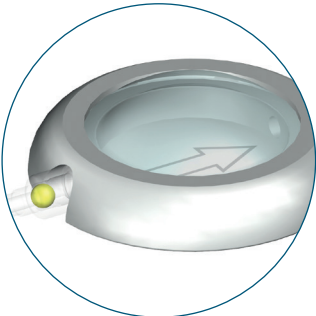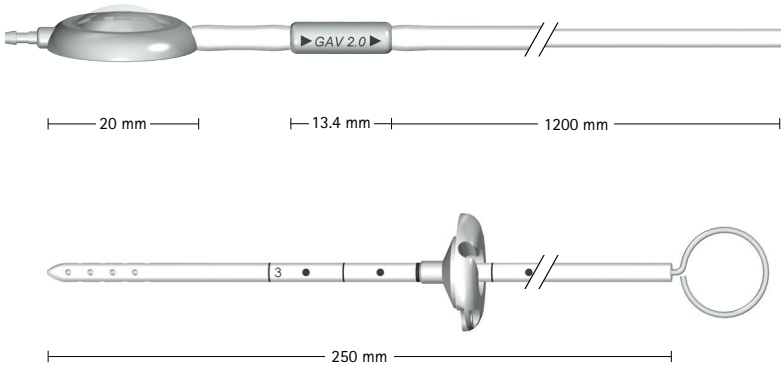

CONTROL RESERVOIR*

Art. No.	Lying	Upright
FX146T	5 cmH ₂ O	20 cmH ₂ O
FX147T	5 cmH ₂ O	25 cmH ₂ O
FX148T	5 cmH ₂ O	30 cmH ₂ O
FX149T	5 cmH ₂ O	35 cmH ₂ O
FX150T	10 cmH ₂ O	25 cmH ₂ O
FX151T	10 cmH ₂ O	30 cmH ₂ O

For pressure level recommendation see page 9.

GAV® 2.0 SHUNT SYSTEM WITH PEDIATRIC CONTROL RESERVOIR

- GAV® 2.0 Valve with integrated pediatric *CONTROL RESERVOIR** and distal catheter (1200 mm)



- Ventricular catheter with pediatric deflector and introducing stylet (250 mm)

* An additional valve in the base of the pediatric *CONTROL RESERVOIR* makes it possible to flush the fluid only in the distal direction. This feature allows for controlling the patency of the ventricular catheter and the distal drainage.

Valve: $d_o = 4.2 \text{ mm}$

Connector: $d_o = 1.9 \text{ mm}$

Catheter: $d_i = 1.2 \text{ mm}$, $d_o = 2.5 \text{ mm}$

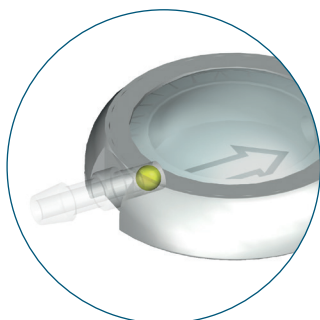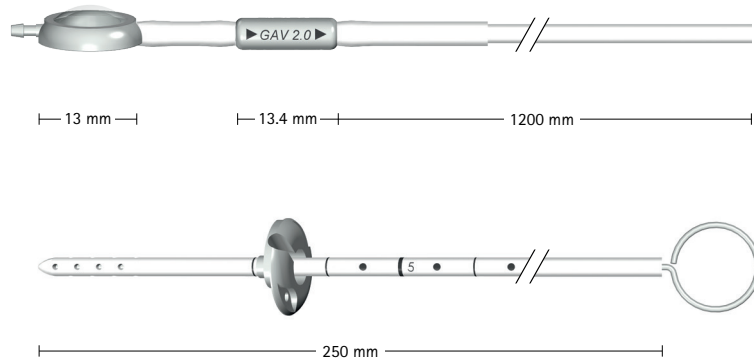

pediatric *CONTROL RESERVOIR**

Art. No.	Lying	Upright
FX152T	5 cmH ₂ O	20 cmH ₂ O
FX153T	5 cmH ₂ O	25 cmH ₂ O
FX154T	5 cmH ₂ O	30 cmH ₂ O
FX155T	5 cmH ₂ O	35 cmH ₂ O
FX156T	10 cmH ₂ O	25 cmH ₂ O
FX157T	10 cmH ₂ O	30 cmH ₂ O

For pressure level recommendation see page 9.

GAV[®] 2.0

GAV[®] 2.0 SHUNT SYSTEM WITH *SPRUNG RESERVOIR*





- GAV[®] 2.0 Valve
with distal catheter (1200 mm)

- *SPRUNG RESERVOIR**
with distal catheter (600 mm)

- Ventricular catheter with
introducing stylet (180 mm)

* An additional valve in the base of the *SPRUNG RESERVOIR* makes it possible to flush the fluid only in the distal direction. This feature allows for controlling the patency of the ventricular catheter and the distal drainage.

Valve: d_o = 4.2 mm
Connector: d_o = 1.9 mm
Catheter: d_i = 1.2 mm, d_o = 2.5 mm

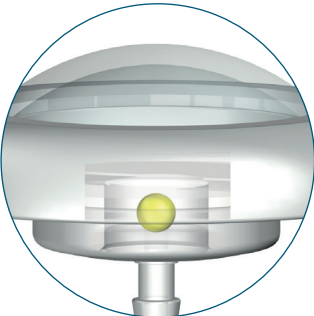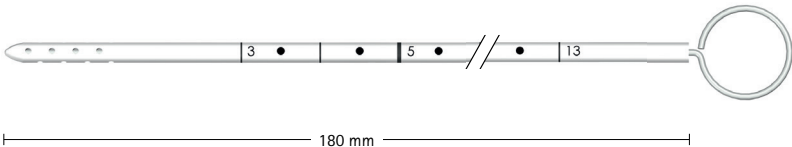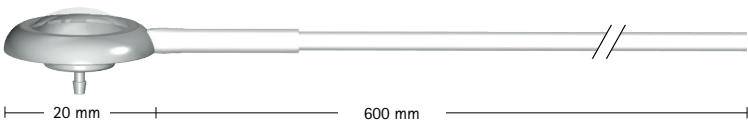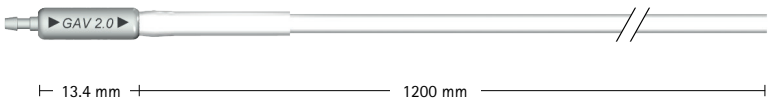


SPRUNG RESERVOIR*

Art. No.	Lying	Upright
FX270T	5 cmH ₂ O	20 cmH ₂ O
FX271T	5 cmH ₂ O	25 cmH ₂ O
FX272T	5 cmH ₂ O	30 cmH ₂ O
FX273T	5 cmH ₂ O	35 cmH ₂ O
FX274T	10 cmH ₂ O	25 cmH ₂ O
FX275T	10 cmH ₂ O	30 cmH ₂ O


For pressure level recommendation see page 9.

GAV® 2.0 SHUNT SYSTEM WITH PEDIATRIC *SPRUNG RESERVOIR*


- GAV® 2.0 Valve
with distal catheter (1200 mm)

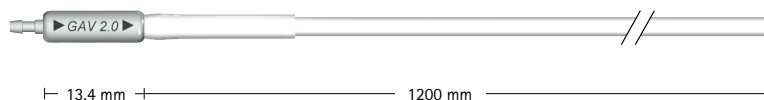

- Pediatric *SPRUNG RESERVOIR**
with distal catheter (600 mm)

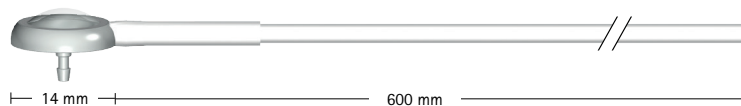

- Ventricular catheter with
introducing stylet (180 mm)

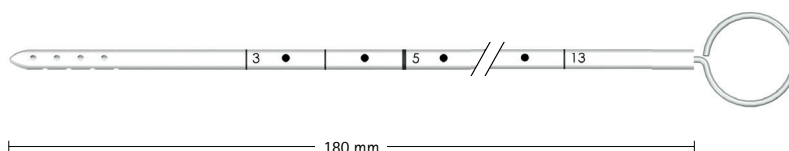


* An additional valve in the base of the pediatric *SPRUNG RESERVOIR* makes it possible to flush the fluid only in the distal direction. This feature allows for controlling the patency of the ventricular catheter and the distal drainage.

Valve: $d_o = 4.2$ mm

Connector: $d_o = 1.9$ mm

Catheter: $d_i = 1.2$ mm, $d_o = 2.5$ mm

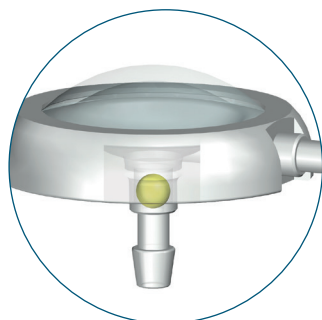

pediatric *SPRUNG RESERVOIR**

Art. No.	Lying	Upright
FX276T	5 cmH ₂ O	20 cmH ₂ O
FX277T	5 cmH ₂ O	25 cmH ₂ O
FX278T	5 cmH ₂ O	30 cmH ₂ O
FX279T	5 cmH ₂ O	35 cmH ₂ O
FX280T	10 cmH ₂ O	25 cmH ₂ O
FX281T	10 cmH ₂ O	30 cmH ₂ O

For pressure level recommendation see page 9.



GAV[®] 2.0

GAV[®] 2.0 SHUNT SYSTEM WITH *SPRUNG RESERVOIR*

- GAV[®] 2.0 Valve with integrated *SPRUNG RESERVOIR** and distal catheter (1200 mm)
- Ventricular catheter with introducing stylet (180 mm)

* An additional valve in the base of the *SPRUNG RESERVOIR* makes it possible to flush the fluid only in the distal direction. This feature allows for controlling the patency of the ventricular catheter and the distal drainage.

Valve: $d_o = 4.2\text{ mm}$
Connector: $d_o = 1.9\text{ mm}$
Catheter: $d_i = 1.2\text{ mm}$, $d_o = 2.5\text{ mm}$

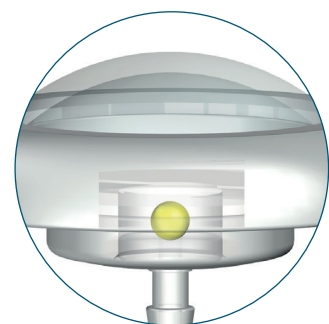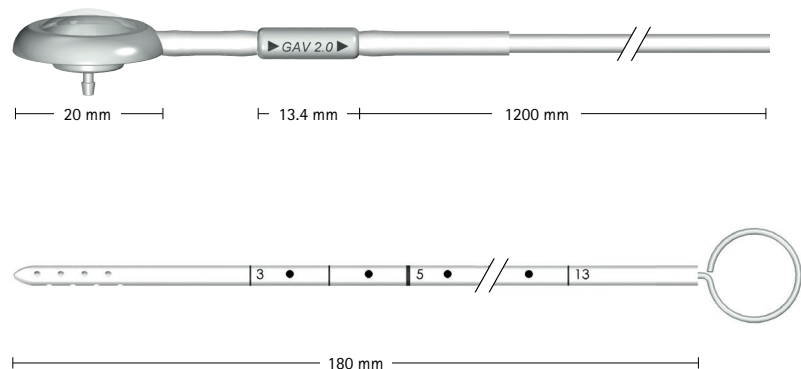

*SPRUNG RESERVOIR**

Art. No.	Lying	Upright
FX170T	5 cmH ₂ O	20 cmH ₂ O
FX171T	5 cmH ₂ O	25 cmH ₂ O
FX172T	5 cmH ₂ O	30 cmH ₂ O
FX173T	5 cmH ₂ O	35 cmH ₂ O
FX174T	10 cmH ₂ O	25 cmH ₂ O
FX175T	10 cmH ₂ O	30 cmH ₂ O


For pressure level recommendation see page 9.

GAV® 2.0 SHUNT SYSTEM WITH PEDIATRIC *SPRUNG RESERVOIR*

- GAV® 2.0 Valve with integrated pediatric *SPRUNG RESERVOIR** and distal catheter (1200 mm)



- Ventricular catheter with introducing stylet (180 mm)

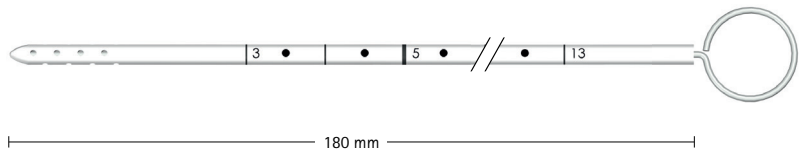


* An additional valve in the base of the pediatric *SPRUNG RESERVOIR* makes it possible to flush the fluid only in the distal direction. This feature allows for controlling the patency of the ventricular catheter and the distal drainage.

Valve: $d_o = 4.2$ mm

Connector: $d_o = 1.9$ mm

Catheter: $d_i = 1.2$ mm, $d_o = 2.5$ mm

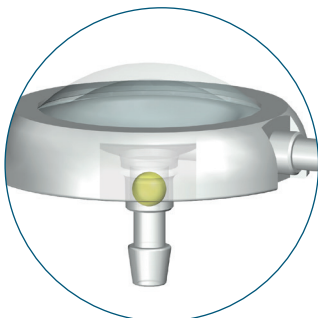

pediatric *SPRUNG RESERVOIR**


Art. No.	Lying	Upright
FX176T	5 cmH ₂ O	20 cmH ₂ O
FX177T	5 cmH ₂ O	25 cmH ₂ O
FX178T	5 cmH ₂ O	30 cmH ₂ O
FX179T	5 cmH ₂ O	35 cmH ₂ O
FX180T	10 cmH ₂ O	25 cmH ₂ O
FX181T	10 cmH ₂ O	30 cmH ₂ O

For pressure level recommendation see page 9.


GAV[®] 2.0

GAV[®] 2.0 SHUNT SYSTEM WITH PEDIATRIC BURRHOLE RESERVOIR


- GAV[®] 2.0 Valve
with distal catheter (1200 mm)

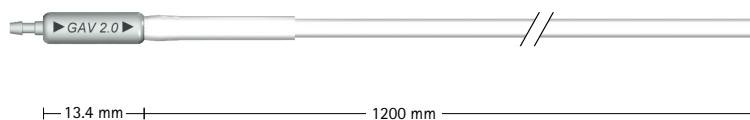

- Pediatric burrhole reservoir
with distal catheter (600 mm)

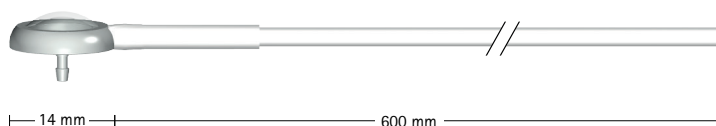

- Ventricular catheter with
introducing stylet (180 mm)

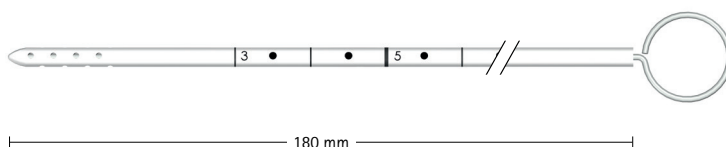


Valve: $d_o = 4.2$ mm
Connector: $d_o = 1.9$ mm
Catheter: $d_i = 1.2$ mm, $d_o = 2.5$ mm

Art. No.	Lying	Upright
FX264T	5 cmH ₂ O	20 cmH ₂ O
FX265T	5 cmH ₂ O	25 cmH ₂ O
FX266T	5 cmH ₂ O	30 cmH ₂ O
FX267T	5 cmH ₂ O	35 cmH ₂ O
FX268T	10 cmH ₂ O	25 cmH ₂ O
FX269T	10 cmH ₂ O	30 cmH ₂ O


For pressure level recommendation see page 9.

GAV® 2.0 SHUNT SYSTEM

- GAV® 2.0 Valve
with distal catheter (1200 mm)

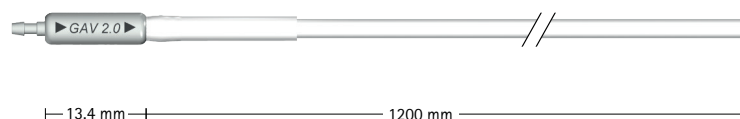

- Ventricular catheter with
pediatric deflector and
introducing stylet (250 mm)

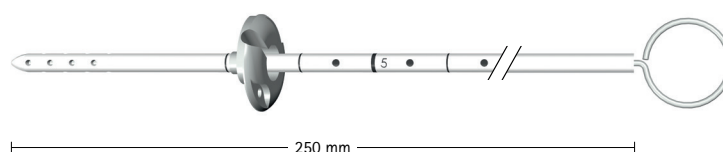

Valve: $d_o = 4.2$ mm
Connector: $d_o = 1.9$ mm
Catheter: $d_i = 1.2$ mm, $d_o = 2.5$ mm


Art. No.	Lying	Upright
FX204T	5 cmH ₂ O	20 cmH ₂ O
FX205T	5 cmH ₂ O	25 cmH ₂ O
FX206T	5 cmH ₂ O	30 cmH ₂ O
FX207T	5 cmH ₂ O	35 cmH ₂ O
FX208T	10 cmH ₂ O	25 cmH ₂ O
FX209T	10 cmH ₂ O	30 cmH ₂ O

For pressure level recommendation see page 9.

GAV[®] 2.0 LP

GAV[®] 2.0 LP, STRAIGHT

- GAV[®] 2.0 LP Valve (straight)
with distal catheter (1200 mm)

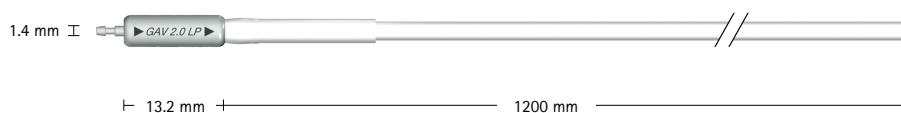


Valve: $d_o = 4.2$ mm
 Connector: $d_o = 1.4$ mm
 for connection with lumbar catheter
 Connector: $d_o = 1.9$ mm
 Catheter: $d_i = 1.2$ mm, $d_o = 2.5$ mm

Art. No.	Lying	Upright
FX222T	5 cmH ₂ O	20 cmH ₂ O
FX223T	5 cmH ₂ O	25 cmH ₂ O
FX224T	5 cmH ₂ O	30 cmH ₂ O
FX225T	5 cmH ₂ O	35 cmH ₂ O
FX226T	10 cmH ₂ O	25 cmH ₂ O
FX227T	10 cmH ₂ O	30 cmH ₂ O

For pressure level recommendation see page 9.

GAV® 2.0 LP, U-FORM

- GAV® 2.0 LP Valve (U-Form)
with distal catheter (1200 mm)

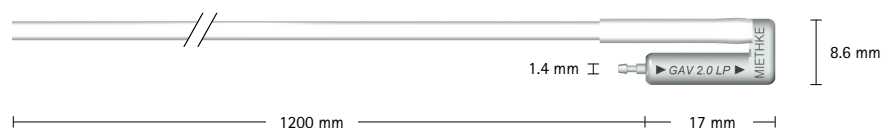


Valve: $d_o = 4.2$ mm
 Connector: $d_o = 1.4$ mm
 for connection with lumbar catheter
 Connector: $d_o = 1.9$ mm
 Catheter: $d_i = 1.2$ mm, $d_o = 2.5$ mm

Art. No.	Lying	Upright
FX228T	5 cmH ₂ O	20 cmH ₂ O
FX229T	5 cmH ₂ O	25 cmH ₂ O
FX230T	5 cmH ₂ O	30 cmH ₂ O
FX231T	5 cmH ₂ O	35 cmH ₂ O
FX232T	10 cmH ₂ O	25 cmH ₂ O
FX233T	10 cmH ₂ O	30 cmH ₂ O

For pressure level recommendation see page 9.

OUR PRODUCTS – YOUR SELECTION

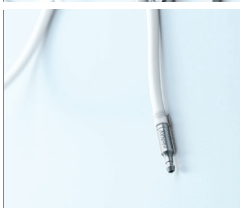	<i>M. blue</i> [®]	<i>proGAV</i> [®] 2.0	<i>GAV</i> [®] 2.0	<i>SHUNT-ASSISTANT</i> [®] 2.0	<i>DUALSWITCH VALVE</i>	<i>miniNAV</i> [®]	Accessories
							
Description	Adjustable gravitational unit with integrated differential pressure valve unit	Adjustable differential pressure valve with gravitational unit	Gravitational valve for treating hydrocephalus	"Add-on" gravitational valve for preventing complications due to excess drainage	Gravitational valve with large flow volumes for CSF	Differential pressure valve, specifically for premature babies and newborns or bedridden, non-mobile patients	
Indication							
LP			✓	✓	✓		
NPH	✓	✓	✓	✓	✓		
Ped. HC	✓	✓	✓	✓		✓	
Adult HC	✓	✓	✓	✓	✓	✓	
Patient							
Bed ridden	✓					✓	
Active	✓	✓	✓	✓	✓	*	
Characteristic							
3-Tesla MR Conditional	✓	✓	✓	✓	✓	✓	
Gravitational unit	✓	✓	✓	✓	✓		
Adjustable	✓	✓					

* in conjunction with *SHUNTASSISTANT*[®] 2.0 or *M.Blue*

Manufacturer acc. to MDD 93/42/EEC

■ **CHRISTOPH MIETHKE GMBH & CO. KG**

Christoph Miethke GmbH & Co. KG | Ulanenweg 2 | 14469 Potsdam | Germany
Phone +49 331 62083-0 | Fax +49 331 62083-40 | www.miethke.com

B. Braun Australia Pty Ltd | Norwest Business Park NSW, 2153 | Tel. 1800 251 705 | info.au@bbraun.com | bbraun.com.au

B. Braun New Zealand Pty Ltd | Level 8, 139 Quay Street, Auckland, 1010 | Tel. 0800 227 286 | info.nz@bbraun.com | bbraun.co.nz

AE 1702 E0821 Rev A 04/26 | ©2026 B. Braun Australia Pty Ltd. ABN 56 002 945 155. All rights reserved. The product trademarks "GAV", "M.blue", "M.blue plus", "miniNAV", "proGAV", "proSA" and "SHUNTASSISTANT" are registered trademarks of Christoph Miethke GmbH & Co. KG in the majority of the world. For more detailed information please contact us: info@miethke.com. Aescuap® is a registered trademark of Aesculap AG. This brochure is intended for use by healthcare professionals only. Subject to technical changes. It is provided for informational purposes and does not constitute medical advice. Always refer to the product's Instructions for Use and consult a qualified healthcare provider for clinical decisions. This material may not be copied, reproduced, or distributed without express written authority of B. Braun Australia Pty Ltd.