




1. Disassembly

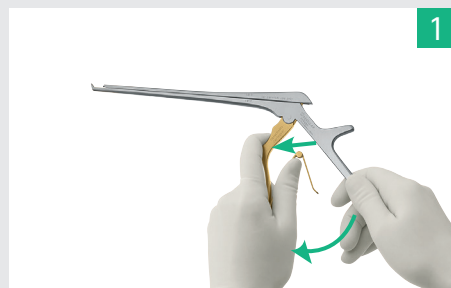

Unhook the spring from the rear handle

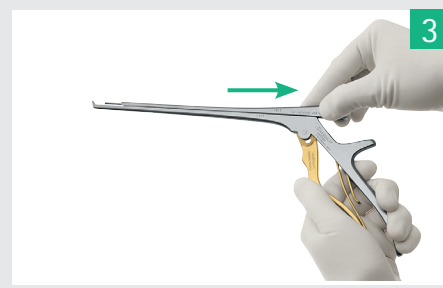

Pull back the slider to remove

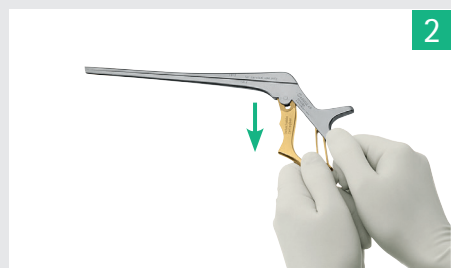

Close the handle completely and pull the inner handle downwards


- Over 80 detachable KERRISON models
- 130° / 90° upbiting
- With standard / thin footplate
- Standard / large jaw punches

2. Cleaning




Please note the instructions for use TA010668

3. Assembly

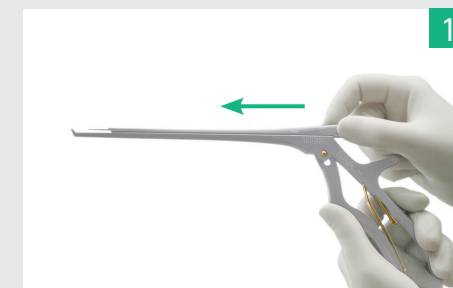

Attach the slider and push it forwards

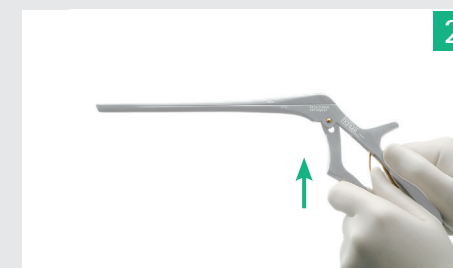

Push the inner handle upwards and open the handle




Insert the spring in the rear handle

- Disassembly and assembly without tools
- Easier verification of cleaning success due to less overlapping surfaces

For further information please see our brochures:


- C26202 Proper Maintenance of Instruments
- C84802 Bone Punches
- C63402 Validated Reprocessing Procedures

4. Maintenance

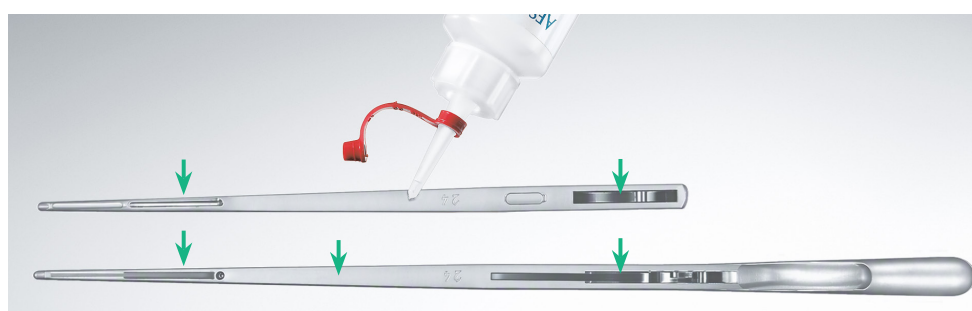


Moving parts (hinges, latches and contact area of sliding part) should be lubricated with a sterilizable and moisture-permeable light weight machine oil such as Sterilit® maintenance oil JG598.

5. Function control

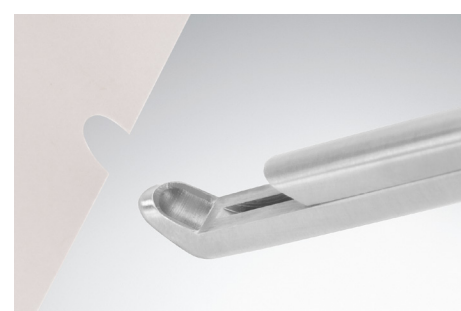

For further information please visit the Aesculap Technical Service

6. Storage

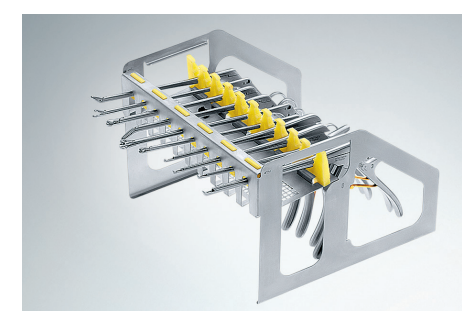


Keep instrument in appropriate storage tray (e.g. storage tray for bone punches JF120R)

7. Sterilization



Please note the instructions for use TA010668

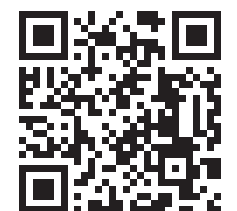

Instruction for Use
<https://eifu.bbraun.com/TA010668>



Aesculap Technical Service
<https://www.bbraun.com/ats>