

AESCULAP®

AdTec® Needle Holder

AESCULAP®

AdTec® Needle Holder







With these needle holders Aesculap® offers a range of instruments specially developed for suturing and knotting. The axial handle allows the instruments to be rotated very easily.

- Specific grip on curved and straight needles; grip on the suture thread.
- An integrated overload protection keeps the holding force of the jaw parts constant. Thus, damage to the needle holder, as well as the needle, can be reduced.
- Simple lock with just two fixed positions following the push ballpoint principle. The needle can thus be gripped with the lock either engaged or released.
- Versatile needle holders available in different diameters and lengths.
- All Durogrip needle holders with tungsten carbide inserts in jaw parts (except PL406R and PL414R).

AESCULAP®

AdTec® Needle Holder




	Description	Diameter	Length	Catalog Number
 1:1	DUROGRIP Needle holder straight	3.5 mm	200 mm	PL419R
		3.5 mm	290 mm	PL418R
		5 mm	310 mm	PL407R
		5 mm	370 mm	PL404R
		5 mm	420 mm	PL416R
 1:1	DUROGRIP Needle holder curved left	5 mm	310 mm	PL408R
		5 mm	370 mm	PL405R
		5 mm	420 mm	PL417R
 1:1	DUROGRIP Needle holder curved right	5 mm	310 mm	PL413R
		5 mm	370 mm	PL415R
 1:1	DUROGRIP Assist. needle holder, 2 x 3 teeth	5 mm	310 mm	PL409R



AESCULAP®

AdTec® Needle Holder

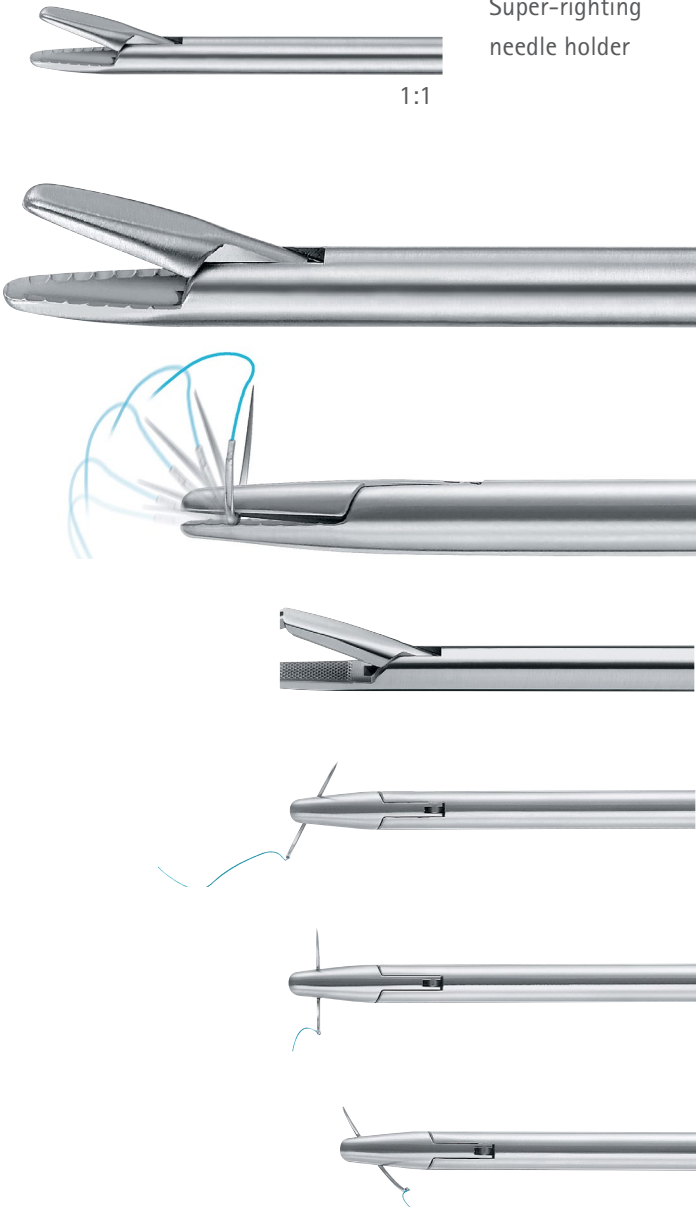


	Description	Diameter	Length	Catalog Number
	Super-righting needle holder	5 mm	310 mm	PL404R

Super-righting needle holder.

Allows automatic upright positioning of the needle at different angles to the instrument shaft.

Additional instrument with same handle

	Description	Diameter	Length	Catalog Number
	Anvil Grasper	10 mm	370 mm	PL406R

AESCULAP® – a B. Braun brand

B. Braun Australia Pty Ltd | Norwest Business Park NSW, 2153 | Tel. 1800 251 705 | info.au@bbraun.com | bbraun.com.au
B. Braun New Zealand | Level 8, 139 Quay Street, Auckland, 1010 | Customer Care 0800 227 286 | Fax (09) 373 5601 | www.bbraun.co.nz

BAUS AU E1591 05/26 ©2026 B. Braun Australia Pty Ltd. ABN 56 002 945 155. All rights reserved. Aesculap® is a registered trademark of B. Braun Melsungen AG. This brochure is intended for use by healthcare professionals only. Subject to technical changes. It is provided for informational purposes and does not constitute medical advice. Always refer to the product's Instructions for Use and consult a qualified healthcare provider for clinical decisions. This material may not be copied, reproduced, or distributed without express written authority of B. Braun Australia Pty Ltd.