



Surecan® Safety II

Power Injectable Safety Needles
for Access Ports

Surecan® Safety II Needles for Access Ports



Surecan® Safety II is the power injectable port access needle with a safety shield that has been designed to offer protection from needlestick injuries, helping to reduce the risk of acquiring a blood borne virus.

Indicated for up to seven days for:

- Infusion of chemotherapy
- Total parenteral nutrition (TPN)
- Antibiotic therapy
- Transfusion
- Pain relief and drainage
- Blood sampling
- Power injections of contrast media

Power Injections

- Suitable for power injections of contrast media up to 325 psi

Coloured clamp

- Possible interruption of infusion at any time
- Easy to identify gauge of the needle

Ergonomic Low Profile Wings

- Provide secure handling during insertion and removal

Extension tubing

- Transparent anti-kink tubing

Special Needle Tip

- Aids needle retention in the septum during power injections

Transparent Base & Non Absorbent Foam Pad

- Helps to reduce risk of infection as fluids are not absorbed into the pad
- Increases patient comfort
- Low profile (1.5 mm) helps to reduce risk of needle movement
- Visualisation of the insertion point

User Safety

- Safety shield is activated when the cannula is extracted from the port.
- An audible click and the green dot indicates that the cannula sharp is locked inside the housing.



Non-Coring Bevel

- Designed to help minimise septum damage

Surecan® Safety II with Y-Site and Caresite® Needlefree Access

Surecan® Safety II is also available with an integrated Caresite needlefree valve which minimises handling steps, helping to reduce the risk of contamination.


Positive Displacement

Upon disconnection of a flush syringe from the Caresite® needlefree valve, an automatic bolus of 0.03 mls of flush solution is delivered which helps to prevent blood reflux and promote catheter patency

LATEX
FREE

MRI
CONDITIONAL

DEHP
FREE

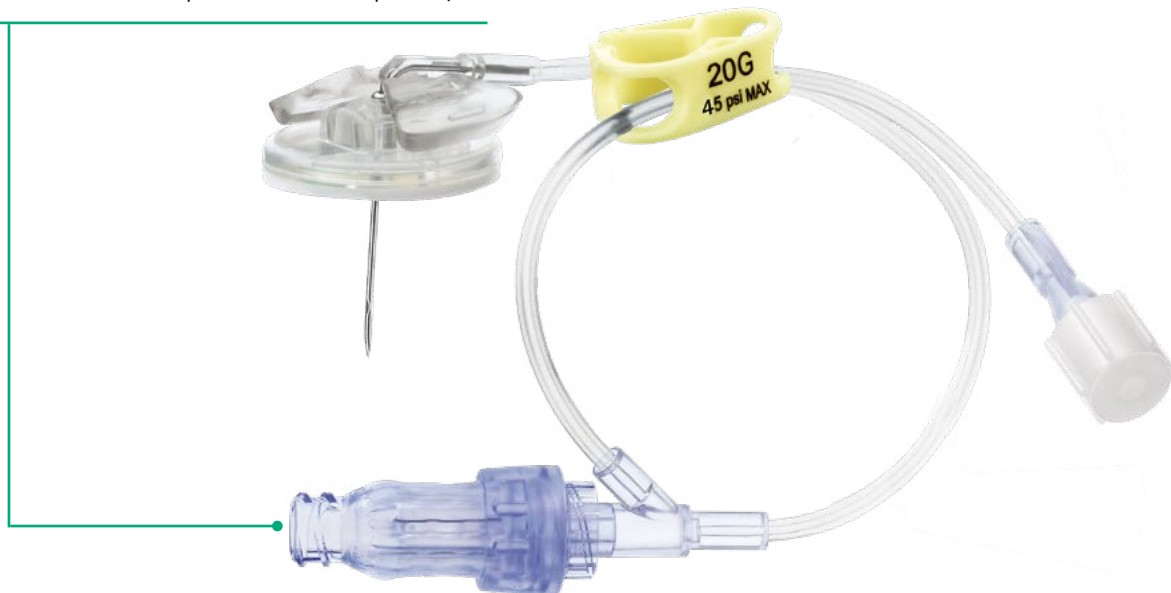

Easy to Clean

The smooth, flat access surface permits thorough and easy disinfection, helping to prevent microbial ingress¹.

Easy to Inspect


The clear outer housing provides complete visibility of the fluid path, allowing for ease of inspection and confirmation of successful flushing.

Easy to Connect

The easy-grip barrel is ergonomically designed to minimise finger slips and the risk of Key-Part touch contamination.

High Pressure (HP)

High Pressure (HP) injection should be administered through the Caresite® needlefree valve on the primary line, which is suitable for use up to 325 psi. Please note that the side arm with the clamp should not be used for HP injection, as this is only suitable for use up to 45 psi.

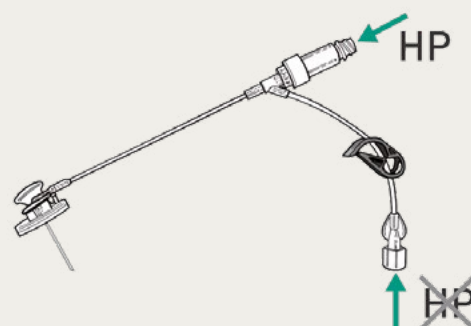

Surecan® Safety II

Size	Length	Priming Volume	Quantity (Per Box)	Reference
19G	12 mm	0.30 ml	20	4447042
19G	15 mm	0.30 ml	20	4447000
19G	20 mm	0.30 ml	20	4447001
19G	25 mm	0.30 ml	20	4447002
19G	32 mm	0.30 ml	20	4447003
19G	38 mm	0.30 ml	20	4447004
20G	12 mm	0.25 ml	20	4447043
20G	15 mm	0.25 ml	20	4447005
20G	20 mm	0.25 ml	20	4447006
20G	25 mm	0.25 ml	20	4447007
20G	32 mm	0.25 ml	20	4447008
20G	38 mm	0.25 ml	20	4447009
22G	12 mm	0.30 ml	20	4447044
22G	15 mm	0.30 ml	20	4447010
22G	20 mm	0.30 ml	20	4447011
22G	25 mm	0.30 ml	20	4447012
22G	32 mm	0.30 ml	20	4447013

Surecan® Safety II with Y-site and Caresite®

Size	Length	Priming Volume	Quantity (Per Box)	Reference
19G	12 mm	0.65 ml	20	4447057
19G	15 mm	0.65 ml	20	4447045
19G	20 mm	0.65 ml	20	4447046
19G	25 mm	0.65 ml	20	4447047
19G	32 mm	0.65 ml	20	4447048
19G	38 mm	0.65 ml	20	4447049
20G	12 mm	0.55 ml	20	4447058
20G	15 mm	0.55 ml	20	4447050
20G	20 mm	0.55 ml	20	4447051
20G	25 mm	0.55 ml	20	4447052
20G	32 mm	0.55 ml	20	4447053
22G	12 mm	0.45 ml	20	4447059
22G	15 mm	0.45 ml	20	4447054
22G	20 mm	0.45 ml	20	4447055
22G	25 mm	0.45 ml	20	4447056

1. Casey A, Karpanen T, Nightingale P, Elliott T. An in vitro comparison of microbial ingress into 8 different needleless IV access devices. J Infus Nurs. 2015 Jan-Feb;38(1):18-25. doi: 10.1097/NAN.0000000000000082. PMID: 25545971.

B. Braun Australia Pty Ltd | Norwest Business Park NSW, 2153 | Tel. 1800 251 705 | info.au@bbraun.com | bbraun.com.au
 B. Braun New Zealand Pty Ltd | Level 8, 139 Quay Street, Auckland, 1010 | Tel. 0800 227 286 | info.nz@bbraun.com | bbraun.co.nz

BAUS AE E1633 A 04/26 ©2026 B. Braun Australia Pty Ltd. ABN 56 002 945 155. All rights reserved. SURECAN® is a registered trademark of B. Braun Melsungen AG. This brochure is intended for use by healthcare professionals only. Subject to technical changes. It is provided for informational purposes and does not constitute medical advice. Always refer to the product's Instructions for Use and consult a qualified healthcare provider for clinical decisions.