





AESCULAP®
AdTec® needle holder



AESCULAP® AdTec® needle holder

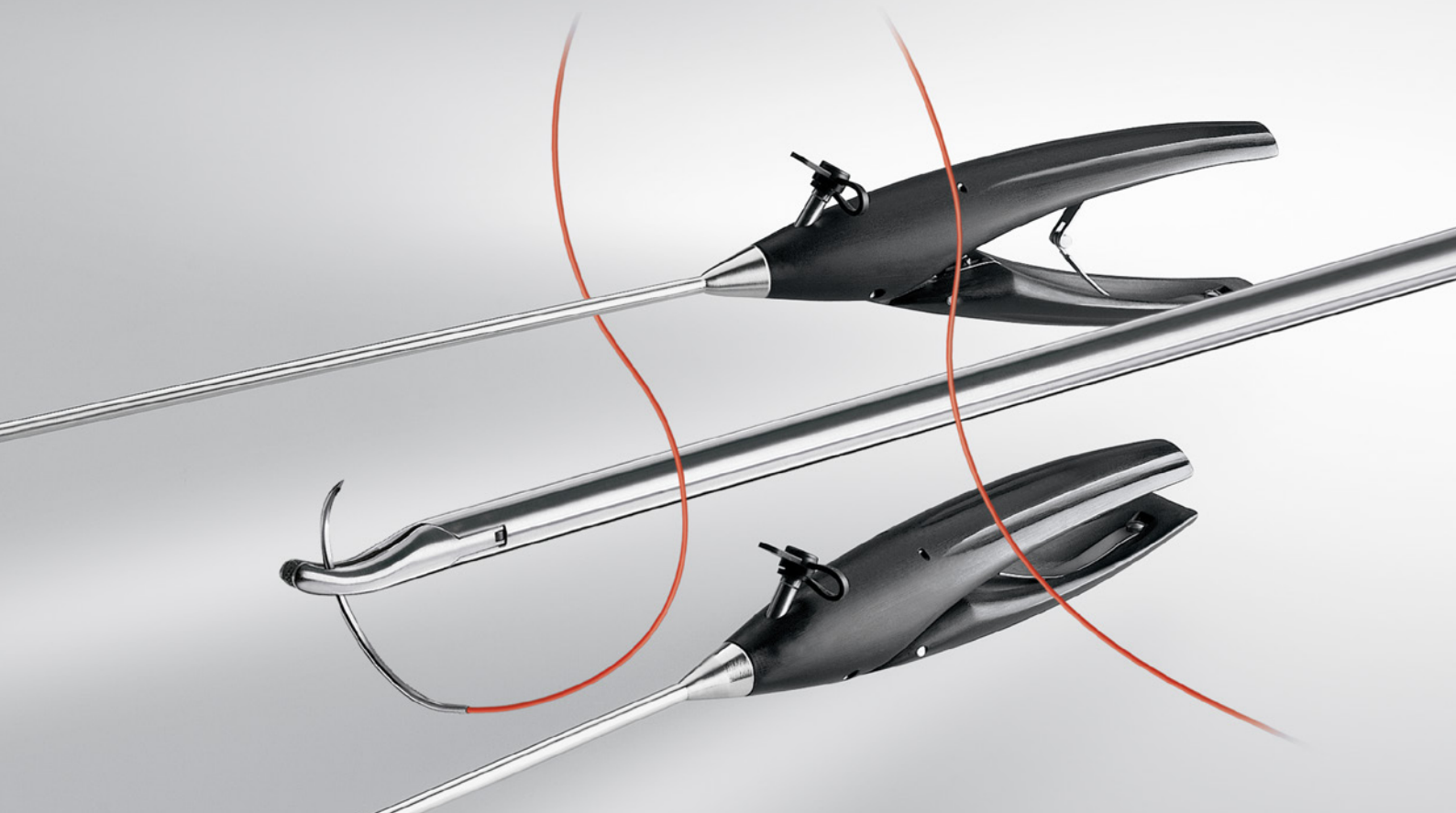

With these needle holders, AESCULAP® offers a range of instruments specially developed for suturing and knotting. The axial handle allows the instruments to be rotated very easily.

- Specific grip on curved and straight needles; grip on the suture thread.
- An integrated overload protection keeps the holding force of the jaw parts constant. Thus, damage to the needle holder, as well as to the needle, can be reduced.
- Simple lock with just two fixed positions, following the push ballpoint principle. The needle can thus be gripped with the lock either engaged or released.
- Improved force transmission to the jaw parts by non-modular instrument design.
- Versatile needle holders available in different diameters and lengths.
- All DUROGRIP needle holders with tungsten carbide inserts in jaw parts (except PL406R and PL414R).

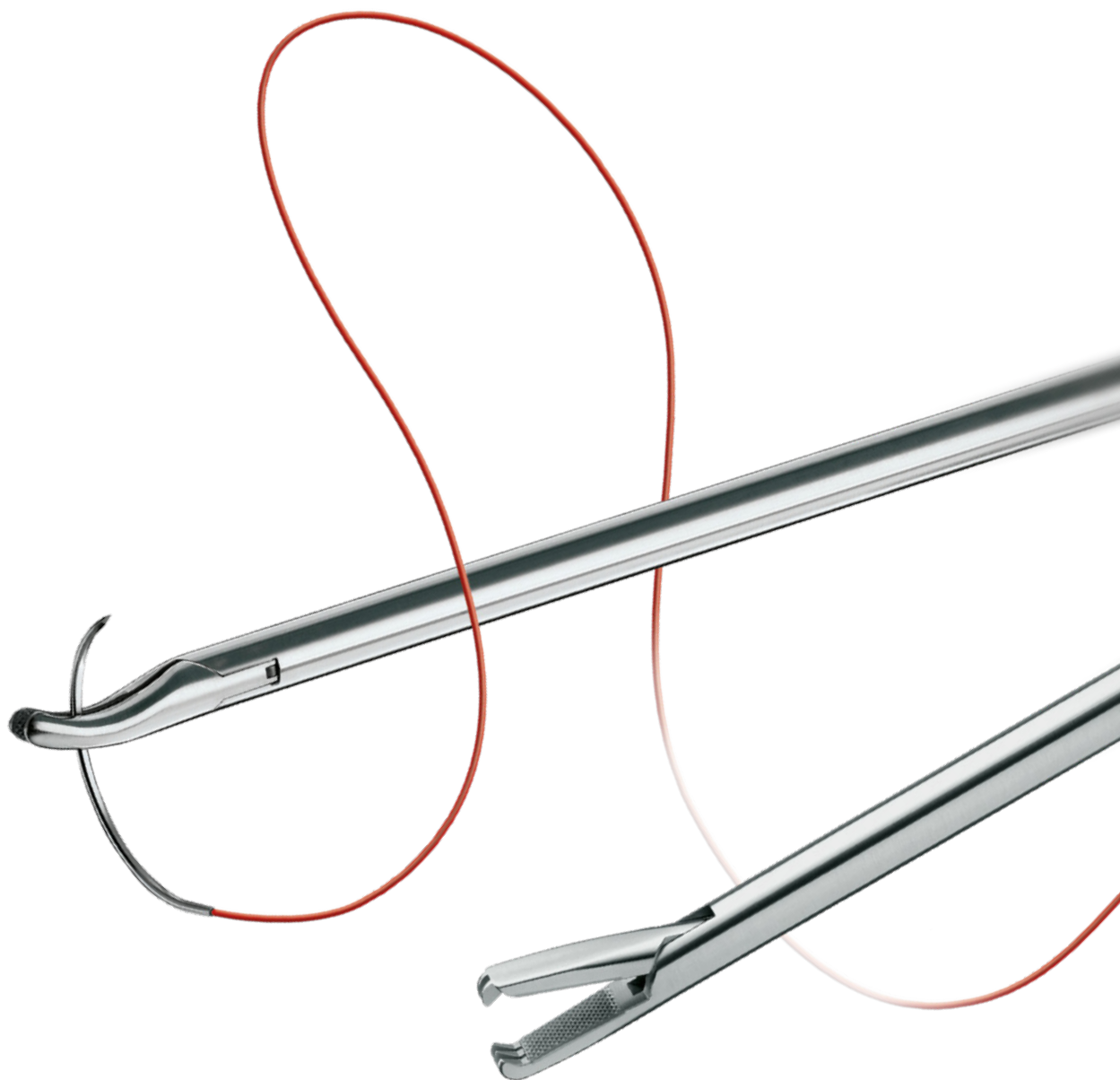


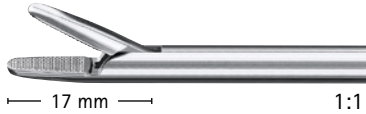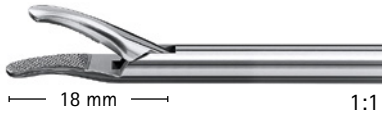Description	Diameter	Length		Catalog number
DUROGRIP	3.5 mm	200 mm		PL419R
Needle holder straight	3.5 mm	290 mm		PL418R
	5 mm	310 mm		PL407R
	5 mm	370 mm		PL404R
	5 mm	420 mm		PL416R
DUROGRIP	5 mm	310 mm		PL408R
Needle holder curved left	5 mm	370 mm		PL405R
	5 mm	420 mm		PL417R
DUROGRIP	5 mm	310 mm		PL413R
Needle holder curved right	5 mm	370 mm		PL415R
DUROGRIP	5 mm	310 mm		PL409R
Assist. needle holder, 2 x 3 teeth				



Description	Diameter	Length		Catalog number
SUPER-righting needle holder	5 mm	310 mm		PL414R

17 mm 1:1

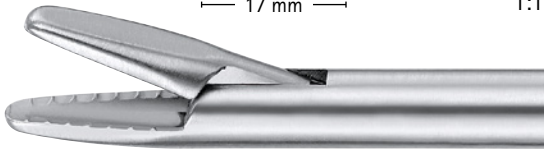

SUPER-righting needle holder.
Allows automatic upright positioning
of the needle at different angles to
the instrument shaft.



Additional instrument with same handle

Description	Diameter	Length		Catalog number
Anvil Grasper	10 mm	370 mm		PL406R

24 mm 1:1

Accessories

Description		Catalog number
Sealing cap irrigation connector, pack à 20 pcs.		EJ751251

AESCULAP® – a B. Braun brand

Aesculap AG | Am Aesculap-Platz | 78532 Tuttlingen | Germany
Phone +49 7461 95-0 | www.aesculap.com

The main product trademark "Aesculap" and the product trademark 'AdTec' are registered trademarks of Aesculap AG.

Subject to technical changes. All rights reserved. This brochure may only be used for the exclusive purpose of obtaining information about our products. Reproduction in any form partial or otherwise is not permitted.