

AESCULAP® Challenger® Ti-P

Pneumatic Reusable Multi-Fire Clip Applier



AESCULAP® Challenger® Ti-P

Pneumatic reusable Multi-fire Clip applicator





- Modular design simplifies disassembly, cleaning, maintenance and sterile preparation
- Innovative technology – pneumatic clip loading
- Reliable vessel occlusion
- Reusable System – reusable applicator and single use clip cartridge
- One System with different shafts

AESCULAP® Challenger® Ti-P

Pneumatic reusable Multi-fire Clip applier



ML

medium-large



for 10 mm clip applier

SM

small-medium



for 5 mm clip applier

Sterile package of clip
cartridge
(clip cartridge + CO₂ cylinder),
last clip color-coded

ML

SM



A diamond-shaped inner clip profile provides an enlarged contact surface




The special closing characteristic of the clip allows a repositioning of the clip and reduces the risk of tissue slippage out of the clip

AESCULAP® Challenger® Ti-P


Pneumatic reusable Multi-fire Clip applier



1:1

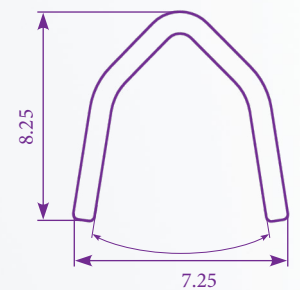

1:1

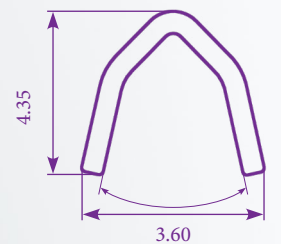

ML

medium-large



Clip Applier	
PL606R	diam. 10 mm, working length 370 mm
PL608R	diam. 10 mm, working length 260 mm
Clip Box	
PL579T	with 12 sterile cartridges of 8 titanium clips and 12 CO ₂ cylinders, last clip color-coded
clip size: medium-large (ML)	



SM

small-medium



Clip Applier	
PL604R	diam. 5 mm, working length 310 mm
PL603R	diam. 5 mm, working length 205 mm
Clip Box	
PL574T	with 12 sterile cartridges of 12 titanium clips and 12 CO ₂ cylinders, last clip color-coded
clip size: small-medium (SM)	



AESCULAP® Challenger® Ti-P


Pneumatic reusable Multi-fire Clip applicator

Components and Spare Parts


ML
medium-large

SM
small-medium

Shaft Schaft	PL536R	PL538R	PL522R	PL609R
complete	shaft length 370 mm	shaft length 260 mm	shaft length 310 mm	shaft length 205 mm
	diam. / Ø 10 mm	diam. / Ø 10 mm	diam. / Ø 5 mm	diam. / Ø 5 mm
Inner Tube Innenrohr	PL536802	PL538802	PL522801	PL601801
	shaft length 370 mm	shaft length 260 mm	shaft length 310 mm	shaft length 205 mm
Push Rod Schubstange	PL536801	PL538801	PL522800	PL601800
	shaft length 370 mm	shaft length 260 mm	shaft length 310 mm	shaft length 205 mm
Handle Griff	PL520R	PL520R	PL520R	PL520R
CO ₂ Cylinder Alone CO ₂ Patrone allein	PL575SU	PL575SU	PL575SU	PL575SU
Pack à 3 pieces, single sterile packed.				
Safety Ring Sicherungsring	PL520301	PL520301	PL520301	PL520301

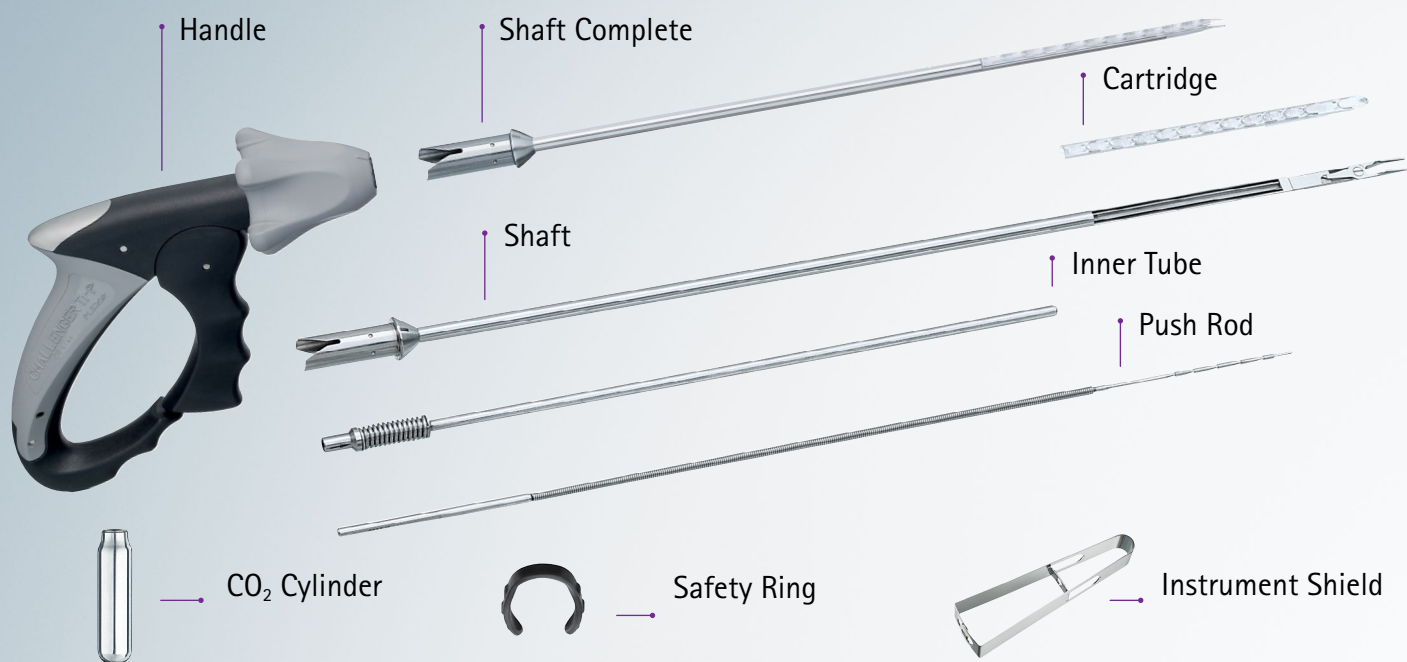

For Care We Recommend

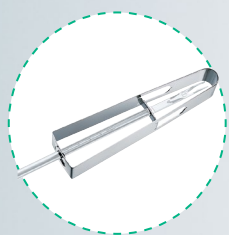


Instrument Shield

PL522810
diam. 5 mm

PL536810
diam. 10 mm

protects jaws during removal, processing and sterile preparation

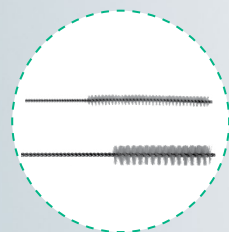

Cleaning Brush

GK469R

for inner tube diam. 5 mm, length 550 mm, pack of 6 pcs.

PM995R

for outer tube diam. 10 mm, length 550 mm, pack of 6 pcs.

Oil Spray

JG600

Sterilit I, Pack of 6 pcs. (not illustrated)

AESCULAP® Challenger® Ti-P

Pneumatic reusable Multi-fire Clip applier

Sterile preparation – In 3 steps ready to use



1. Assembly of handle and shaft



2. Insert CO2 cylinder



3. Insert clip cartridge

AESCULAP® – a B. Braun brand

B. Braun Australia Pty Ltd | Norwest Business Park NSW, 2153 | Tel. 1800 251 705 | info.au@bbraun.com | bbraun.com.au
B. Braun New Zealand | Level 8, 139 Quay Street, Auckland, 1010 | Customer Care 0800 227 286 | Fax (09) 373 5601 | www.bbraun.co.nz

BAUS AU E1594 05/26 ©2026 B. Braun Australia Pty Ltd. ABN 56 002 945 155. All rights reserved. Aesculap® is a registered trademark of B. Braun Melsungen AG. This brochure is intended for use by healthcare professionals only. Subject to technical changes. It is provided for informational purposes and does not constitute medical advice. Always refer to the product's Instructions for Use and consult a qualified healthcare provider for clinical decisions. This material may not be copied, reproduced, or distributed without express written authority of B. Braun Australia Pty Ltd.