
01907/2020/09/036

NaCl 0,9% B. Braun

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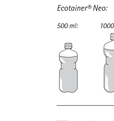en NaCl 0,9% B. Braun Instruction for use for patients or caregivers

Instruction for use for professional users: Please refer to B. Braun homepage: https://eifu.bbraun.com/. Or [QR code]

1000 ml solution contain:
Sodium chloride in Water for Injection: 9.00 g

Electrolytes: mmol/l
Na+: 154
Cl-: 154

Theoretical Osmolarity: 308 mOsm/l
This package unit of sale contains multiple containers of irrigation solution. Such containers shall not be distributed individually unless an instruction for use is provided with each container.

Ecotainer® Neo:
– Irrigation bottle with screw-cap.
– NaCl 0.9% B. Braun is stored in the bottle and can be used without additional devices.

Ecolav® Neo:
– Irrigation bottle with twist-off cap.
– NaCl 0.9% B. Braun is stored in the bottle and can be used without additional devices.

NaCl 0,9% B. Braun:
– Provides mechanical cleansing for sterile irrigation of eyes, body cavities, tissues, wounds, and burns, including during and after surgery.
– Helps maintain a clear surgical field by removing blood and tissue.
– Is used for rinsing medical instruments and equipment as well as for moistening airways by inhalation.
– Electrolyte concentration is similar to plasma.
– Do not reuse the product.

Sterility
Sterile and endotoxin-free irrigation solution. Sterilization is performed by moist heat.
Only use NaCl 0.9% B. Braun if the tamper evidence seal/container is undamaged and the solution is clear.

Intended purpose
NaCl 0.9% B. Braun is a ready for use, sterile, single use medical device, intended to be used for irrigation applications (irrigation solution).

Indication
NaCl 0.9% B. Braun is used as irrigation solutions or moistening agent in patients in need of:
– surgical interventions where electrolyte solutions can be safely used (intra- and postoperative irrigation)
 – rinsing of urinary catheters
 – endoscopic examinations of a body cavity
 – arthroscopy
 – irrigation of wounds and burns
 – moistening wound tamponades, cloths and dressings
 – stoma care
 – nasal irrigation and gargling
 – mechanical irrigation of the eyes
 – moistening the airways by means of inhalation

Patient population
NaCl 0.9% B. Braun can be used for all patients under consideration of the respective indications, contraindications and warnings for use.

Intended User
The products should only be used by trained professionals such as healthcare workers, licensed physicians, nurses, pharmacists, and caregivers.
After proper medical assessment, instruction and documentation by the healthcare worker, patients and caregivers may also handle the products according to national guidelines.

Contra-indications
Do not use NaCl 0.9% B. Braun:
– To inject or infuse
– If you are allergic to any of the substances or materials employed.

Side effects
– Using NaCl 0.9% B. Braun can cause burning sensations, irritations, or discomfort at the application site.
– Displaced catheters or drainage tubes can lead to unintended irrigation or infiltration of body structures.
– Using too much volume or pressure can damage tissues or cause leaks.
– Contamination due to improper handling may cause infection, which can lead to sepsis.
– Continuous or excessive irrigation can cause fluid overload, leading to symptoms like electrolyte imbalances, heart issues, swelling, difficulty breathing, nausea, and kidney problems.

Precautions
– Do not drop the container. If dropped, discard it.
– If you are unsure about using the irrigation solution, consult your general practitioner or pharmacist.
– NaCl 0.9% B. Braun is compatible with adrenaline and noradrenaline and can be admixed directly prior to use.
– In hospital settings only, NaCl 0.9% B. Braun may be warmed up to not more than 40°C for maximum two weeks.

Warning
– During irrigation in surgery or post-operative settings, healthcare practitioners will monitor for excessive fluid absorption. If fluid overload is suspected, irrigation will be stopped, and the patient will be evaluated. Symptoms may appear later. For irrigation of wounds and burns, there is no risk of absorbing the solution into the body if used according to intended use.
– In surgery or post-operative settings, when there is a risk of fluid absorption, NaCl 0.9% must be used with caution in patients with low level of potassium or too high blood levels of sodium and chloride or with conditions that may be associated with it (e.g. impaired kidney function, heart failure, swelling, water on your lungs, high blood pressure). Imbalances should be corrected before the use of NaCl 0,9%.
– NaCl 0.9% must not be used with types of electrosurgery that are not safe with electrolyte solutions.
– Sterile precautions are to be taken when handling the irrigation solution.
– Do not save opened containers; discard any leftover solution.
– Do not use NaCl 0,9% B. Braun after the expiry date on the label.
– Re-using single-use devices can lead to contamination and functional issues, posing risks to patients and users.
– Keep out of sight and reach of children.

Operating Instructions
Ecotainer® Neo:
1. Grip container and open container by removing the twist-off cap.
2. Pour solution in a sterile kidney bowl or address patient area by applying slight pressure on the soft bottle.

Ecolav® Neo:
1. Grip container and open container by removing the twist-off cap.
2. Pour solution in a sterile kidney bowl or address patient area by applying slight pressure on the soft bottle.

Duration of use
– The amount and duration of use depend on the intervention and postoperative recovery. The container does not need to be completely emptied.
– To avoid significant absorption during surgery, your healthcare professional must avoid prolonged or continuous use of the solution.
– In hospital post-operative settings, irrigation may last more than 60 minutes, but not more than 24 hours continuously.

Disposal
Ask your pharmacist how to dispose of the solution and containers no longer required.

Storage and handling conditions
Please refer to the storage information on the packaging.
In hospitals, compatible disinfectants can be used on the surface of irrigation containers.

Clinical benefits B Performances
The clinical benefit of the irrigation solution is related to the result of the associated therapy, when used as intended.

Notice to the user
All serious incidents related to this product must be reported to B. Braun and the competent authority in the country where the product is being used.
Instructions For Use for professional users are available on the B. Braun homepage: https://eifu.bbraun.com/.

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B. Braun Melsungen AG,